

Laguna



Care Relax \

Geriatric reclining armchair finished to an extremely high standard (a step up from the Echo in terms of comfort) offering excellent comfort. The recliner has an ergonomically shaped back, which makes it ideal for people who have to sit down for extended periods of time (elderly who are not mobile and need plenty of care). The Laguna has numerous adjustment/reclining options (which is a plus compared to the Echo). This means that fellow residents or nurses can help put the armchair in the desired position if the user is no longer capable of doing so.

What makes the Laguna Care Relax so unique? \

- \\ Multiple reclining/adjustment mechanisms, also electrically adjustable
- \\ Ergonomically shaped back
- \\ Unique design
- \\ Different side panels
- \\ Comfortable sitting over long periods
- \\ Optional removable headboard
- \\ Maintenance-friendly removable seat

OVERVIEW LAGUNA CARE RELAX



Armrests

Available with 7 different armrests, both open and closed and in different types of materials.

More information on page 5

Adjustment mechanism

5 different adjustment mechanisms are possible, both manual and electric. Available with or without leg support and optional (adjustable) footrest for Laguna without leg support.

More information on pages 3 & 4

Frame & castors

Available with 4 different types of base, depending on the frame, castor size and braking method.


More information on page 5

Options & accessories

Laguna Care Relax has a whole host of options and accessories to adapt the Care Relax to the user's needs.

More information on page 6

RECLINING MECHANISMS

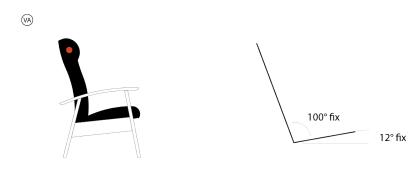


FIXED VERSION


Fixed version, not adjustable, without leg rest. Optional: fixed or adjustable footstool.




sc

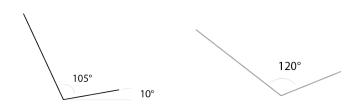


The seat and back can be adjusted steplessly by means of a gas spring, includes posture correction feature. Without leg rest. Optional: fixed or adjustable footstool.

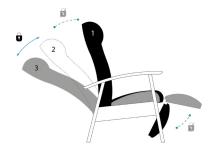

QT

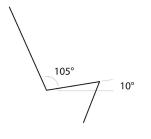

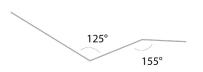


105°

155°


The seat and back can be adjusted and the back can be adjusted separately (by means of a gas spring), includes posture correction feature.

Subject to change without notice, latest update 2021-01, images shown are for illustration purposes only

VERSTELMECHANIEKEN

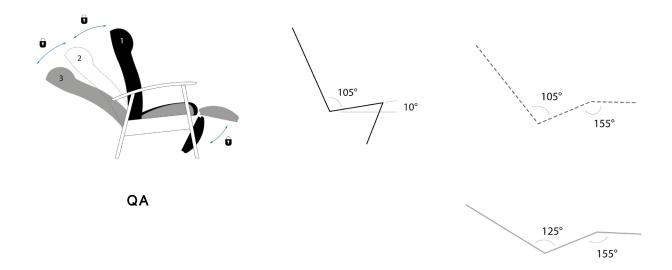

The seat and back can be adjusted and the back can be adjusted separately (by means of a gas spring), includes posture correction feature. Separate foldable leg rest with safety locking system.

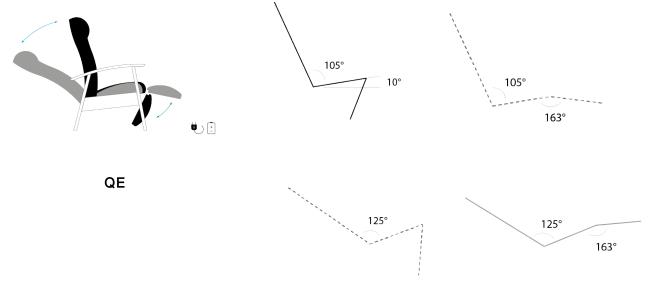

The seat and leg rest can be adjusted electrically (2 motors), includes posture correction feature. Leg rest can be adjusted separately.

Choice: with fixed floor adapter or fixed rechargeable battery.


ARMRESTS

The Laguna can be equipped with various side panels/armrests: from a comfortable, fully padded armrest to a lowerable armrest. Depending on the type, it may also be possible to attach a dining tray to the armrest.



FRAME & CASTORS

The following versions are possible:

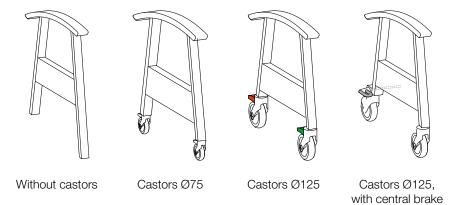

5

TECHNICAL FEATURES

DIMENSIONS

TYPE	В	D	н						kg		*
Laguna	70-82	85-92	115	55	50	50	75	69-70	28-50	yes	***
Footstool fixed	53	43	43	_	_	_	_	_	4,1	yes	***
Footstool adjustable	53	43	43	-	_	_	-	-	4,2	yes	***

 \star = domestic use / \star \star = project-based use / \star \star = intensive project-based use

BACK, SEAT AND LEG REST

Back, seat and legrest (except for the VA and SC, which have no leg rest) are upholstered. Excellent seating comfort. Ergonomically shaped back consisting of a plywood shell and foam padding. The seat is made of plywood with webbing and HR foam. Standard: with removable seat. Optional: with integrated leg rest (but without a removable seat).

OPTIONS & ACCESSORIES

There are numerous options and accessories available for the Laguna, including a dining tray, IV pole... Please have a look at our technical data sheets and product brochure.

Subject to change without notice, latest update 2021-01, images shown are for illustration purposes only

OVERVIEW OF MODELS



LAGUNA CLASSIC FIXED



LAGUNA CLASSIC CASTORS Ø 75



LAGUNA CLASSIC CASTORS Ø 125



LAGUNA HARMONY FIXED



LAGUNA HARMONY CASTORS Ø 75



LAGUNA HARMONY CASTORS Ø 125



LAGUNA WOOD FIXED



LAGUNA WOOD CASTORS Ø 75



LAGUNA WOOD CASTORS Ø 125

OVERVIEW OF MODELS



LAGUNA COSY FIXED



LAGUNA COSY CASTORS Ø 75



LAGUNA COSY CASTORS Ø 125



LAGUNA COSY CASTORS Ø 125 CENTRAL BRAKE



LAGUNA FOREST FIXED



LAGUNA FOREST CASTORS Ø 75



LAGUNA FOREST CASTORS Ø 125



LAGUNA FOREST CASTORS Ø 125 CENTRAL BRAKE



LAGUNA FOREST NS CASTORS Ø 125



LAGUNA FOREST NS CASTORS Ø 125 CENTRAL BRAKE



LAGUNA HOSPITAL NS CASTORS Ø 125



LAGUNA HOSPITAL NS CASTORS Ø 125 CENTRAL BRAKE